



Hanson M1 & M2 Series Workstations

The newest member of the HLS family

HLS M-series bench systems are highly configurable, and offer numerous accessory options designed to meet a variety of project budgets.

The backbone of the M1 & M2 bench is a robust 2" x 3" 14-gauge powder-coated steel upright for maximum stability and extended life.

Sizing choices:

- Widths - 4', 5', 6'
- Depths - 30", 36", 60" (double-sided)
- Height adjustable in 1" increments from 30-3/4 to 36-3/4".

Power, data, and gas service options can be pre-installed into the frame, resulting in a clean, uncompromised work center.

Additional accessories:



- 1" thick lab-grade phenolic work surfaces
- Shelving - multiple depth choices with various retaining lip styles
- Mobile and suspended cabinets - multiple door and drawer configurations
- Bin box mount
- Horizontal power strip



M1 bench (single-sided)



M2 bench (double-sided)

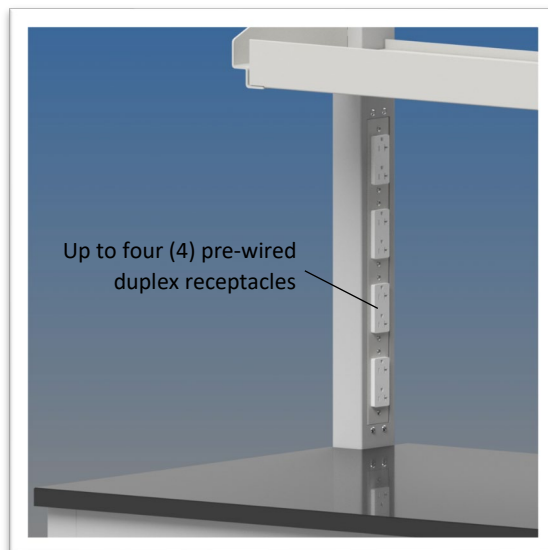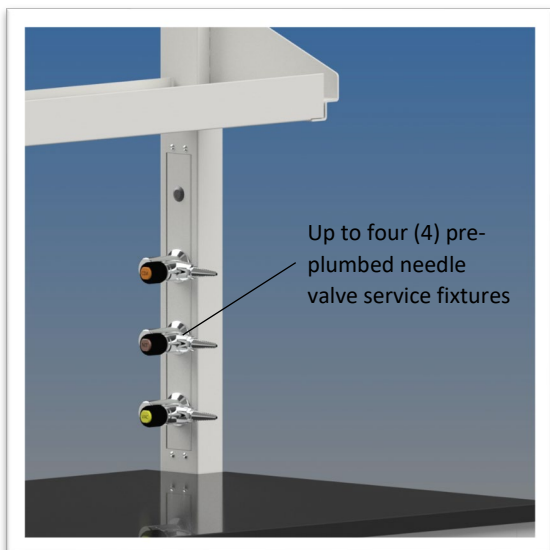




Designed and manufactured in Camarillo, California 

Hanson Lab Solutions - 805.498.3121


www.HANSONLAB.com

HLS M1 & M2 WORKSTATIONS


- **FRAME CONFIGURATIONS**
 - M1: SINGLE-SIDED, 30 OR 36" DEEP
 - M2: DOUBLE-SIDED, 60" DEEP
- **14GA STEEL 2 x 3" REAR POST, 2 x 2" FRONT LEG CONSTRUCTION**
- **SIX-POSITION ADJUSTABLE WORKSURFACE HEIGHT (1" INCREMENTS)**
- **STATIC LOAD CAPACITY, EQUALLY DISTRIBUTED**
 - COMBINED SYSTEM, NO CASTERS : 1200 LBS (880 LBS W/ CASTERS)
 - LOWER SHELF : 200 LBS
 - UPPER SHELVES : 150 LBS
- **AVAILABLE OPTIONS:**
 - UPPER SHELF KIT (SHOWN) IN VARIOUS SIZES / STYLES
 - LOWER SHELF KIT (SHOWN)
 - POWER, DATA AND GAS SERVICE INTEGRATION INTO UPRIGHT POSTS (SHOWN)
 - CHEM RESISTANT PHENOLIC WORK SURFACE (SHOWN)
 - SIDE APRON OPTION TO ACCOMMODATE ALTERNATE DEPTH WORK SURFACES
 - SUSPENDED CABINET (SHOWN) IN VARIOUS SIZES / STYLES
 - MOBILE CABINET IN VARIOUS SIZES / STYLES
 - HORIZONTAL POWER STRIP KIT (SHOWN)
 - BIN BOX MOUNT BRACKET (SHOWN)
 - SWIVEL CASTER KITS (W/BRAKE OR LEVELING OPTION)
- **SEE PRICE LIST FOR FRAME AND ACCESSORY PART NUMBERS**
- **ASSEMBLY REQUIRED**

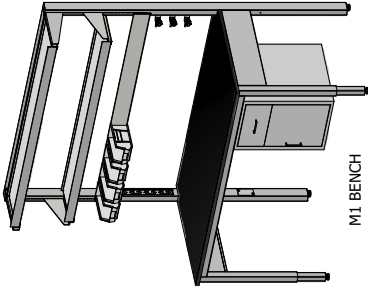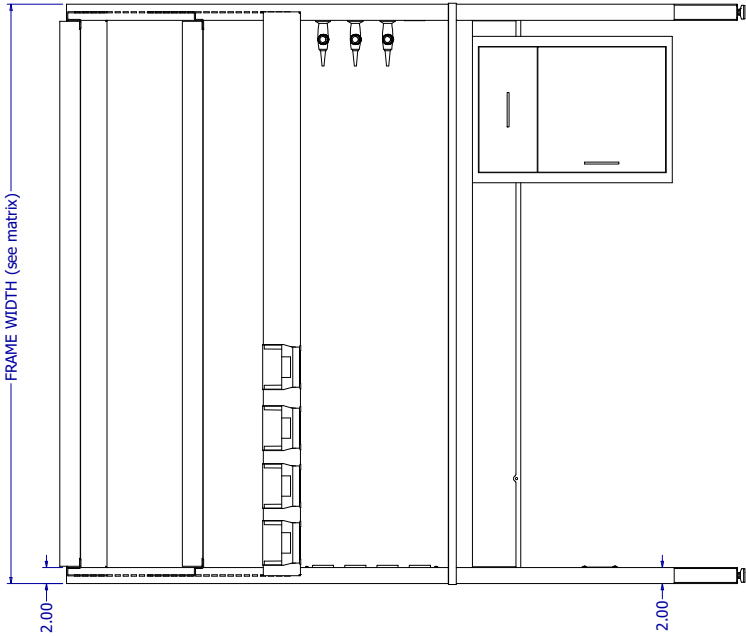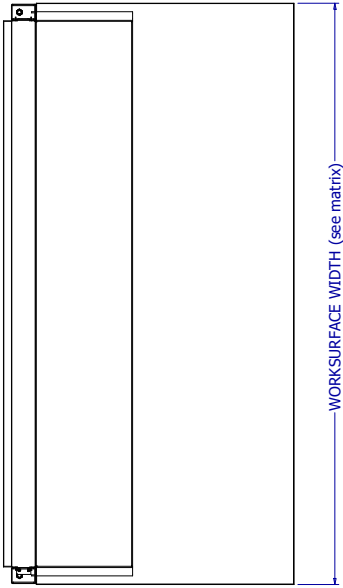

M1 BENCH

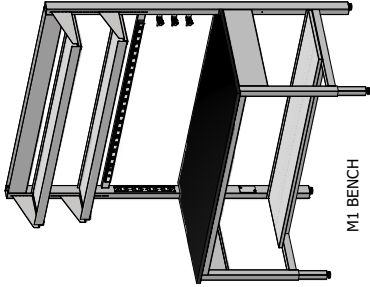

M1 BENCH

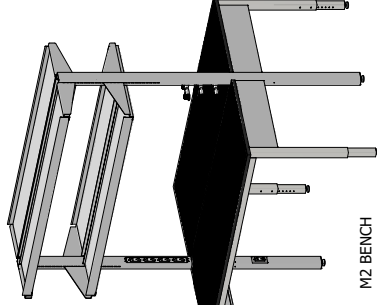

M2 BENCH

SIZE MATRIX		
Nominal (width x depth)	Frame width	Worksurface width
36 x 30, 36 or 60*	35-3/4"	36"
48 x 30, 36 or 60*	47-3/4"	48"
60 x 30, 36 or 60*	59-3/4"	60"
72 x 30, 36 or 60*	71-3/4"	72"

*58" or 60" depth indicates double-sided M2 bench

DATE	4/24/2023	HLS HANSON LAB SOLUTIONS	
DESIGNED	J. HANSON	747 Cole Blvd, Corona, CA 92621 951-483-3121	
QA		TITLE	BENCH ASSY, HLS M1/M2 SERIES 2x3, ADJ HT
APPROVED		FINISH	POWDERCOAT PER JOB REQUIREMENT
UNLESS OTHERWISE NOTED, ALL DIMENSIONS ARE IN INCHES		SIZE	D
		JOB #	
		QUANTITY	
		SCALE	1:1
		PART NO.	SEE PRODUCT LIST
		REV	-
		SHEET 1 OF 1	